

THORATEX 05

Semirigid thoracic support

Ref. SRT 142

Indications

It is used to correct vicious positions of the thoracic spine, or in case of kyphosis.

Description


The support is made of soft textile, which does not harm while wearing, provided with:

- one posterior component that supports the thoracic spine, made of rigid fabric (which does not stretch)
- 2 rigid straps around the shoulders, correcting the body posture by traction. The 2 straps have one elastic and one rigid component, are positioned over the head of the humerus and are adjustable by cutting. To shorten the straps, remove the "crocodile" type fastening from the extremities, cut the strap to the optimum length on the prefigured transverse marks and reposition the hook and loop.
- 2 paravertebral integrated metal rods.


Hook and loop fastener on the front side.

Composition


Two and three-layer laminated textile. The laminated textile ensemble behaves like a moisture regulator, ensuring a quick transfer of the moisture to outward.

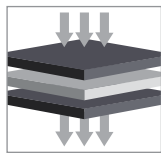

air transfer

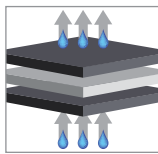

water transfer

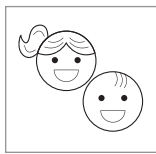

also available in
Junior version




Product details



- hook and loop fastening on the front side



- 2 straps positioned over the head of the humerus for an efficient traction

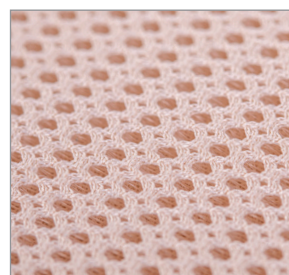

- spatialized 3D mesh



- Available color:
beige

Available sizes

Sizes	1	2	3
Torso circumference (cm)	75-95	95-115	115-135
Reference	SRT142M1	SRT142M2	SRT142M3

The circumference of the torso is measured

