

## MALLEOFIX PLUS

Mobile ankle support with elastic strap

Ref. SRT 312

### Indications

It is used in cricks, luxations, ankle articulation instability.

### Description

The orthosis is made of textile, provided with:

- 2 lateral pads that fit the malleolus (polyethylene foam)
- 2 elastic straps crossed around the ankle for compression and fastening.

The heel area is open.


Hook and loop fastener on the front side.

Height: 21 cm.


The orthosis is an ankle support, and the two elastic straps around the ankle exert additional compression, ensuring stability while wearing.

### Composition


Three-layer laminated textile. The laminated textile ensemble behaves like a moisture regulator, ensuring a quick transfer of the moisture to outward.

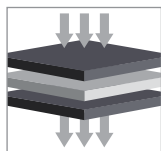

air transfer

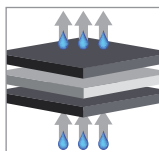

water transfer

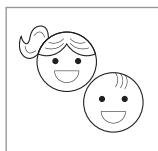

also available in  
Junior version



### Product details



- hook and loop fastener on the front side
- elastic straps crossed around the ankle for compression and fastening



- open heel area




- lateral pads made of polyethylene foam



- Available color:  
black

## Available sizes

Sizes	1	2	3	4	5
<b>Ankle circumference (cm)</b>	20-22	22-24	24-26	26-28	28-30
<b>Reference</b>	<b>SRT312M1</b>	<b>SRT312M2</b>	<b>SRT312M3</b>	<b>SRT312M4</b>	<b>SRT312M5</b>

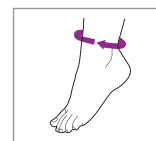

The circumference of the ankle is measured around the malleolus