

TRIACARP JUNIOR

Fixed wrist hand brace

Code: **SRT 222**

Indications

It is used in sprains, dislocations, carpal tunnel syndrome, radial nerve palsy and fractures at the level of the area, both as orthopedic and postoperative treatment.

Description

It is made of textile, provided with:

- an anatomically shaped, mouldable metal palm splint
- extra strip around the wrist for a better fit.

Hook and loop fastener on the dorsal side of the hand.

Length varies between 14-18 cm, depending on the size.

The orthosis stabilizes and immobilizes the radiocarpal joint.

Composition

The orthosis is made of a three-layer laminated textile:


- **Asterix** - The polyamide fiber has a high resistance to chemical and biological agents, resists abrasion and does not absorb water.
- **Moltopren (MTP)** - air-permeable open-cell polyurethane foam (airflow in both directions is inversely proportional to foam density) with high moisture absorbing properties.
- **Jersey cotone SOFT** - 100% cotton or cotton with spandex - elastane (for increased elasticity).

SOFT


The code for the balsamic treatment of the cotton fabric to provide softness to the touch and to avoid skin irritation.

Properties of the fabric

The laminated ensemble behaves like a moisture regulator: the cotton fabric absorbs body sweat and transfers it to the PU foam and outward.

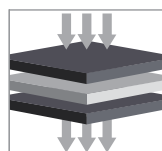

airflow

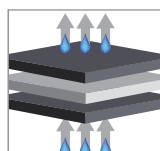

water transfer



also available in adult version



Product details




- strip around the wrist for a better fit



- mouldable metal palm splint

Available sizes, left/right versions

Sizes	1	2	3
Wrist circumference (cm)	10-12	12-14	14-16
Length (cm)	14	16	18
Code	SRT222M1	SRT222M2	SRT222M3

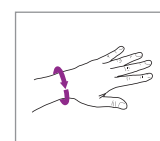

After the code must be added **L for Left, R for Right**

The circumference of the wrist is measured